
GEBRAUCHSANWEISUNG
OPERATING MANUAL
MODE D'EMPLOI
INSTRUCCIONES PARA EL USO
MANUALE D'USO

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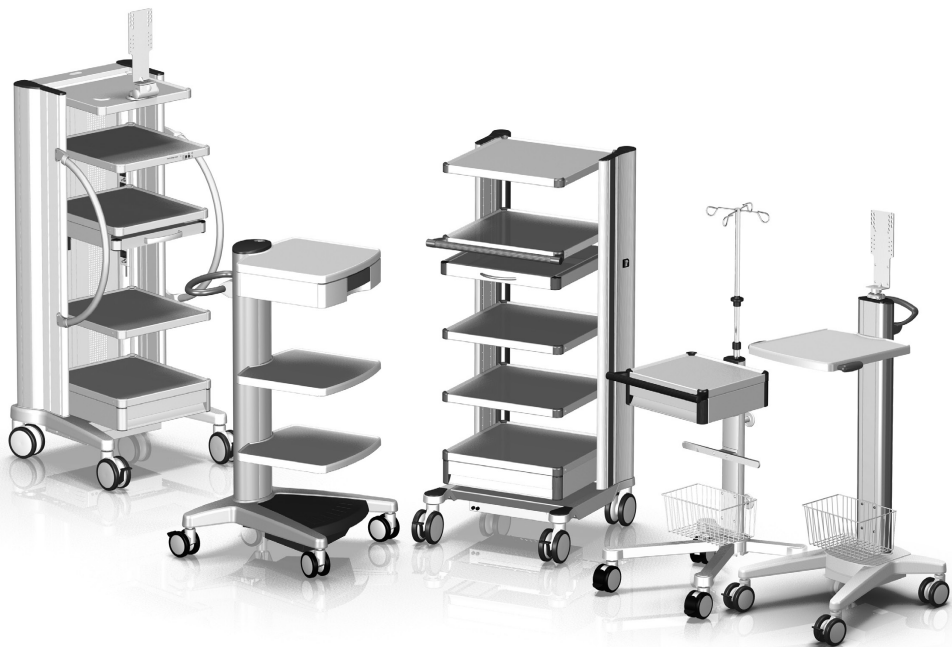
This is a Class I medical product in accordance with the European Medical Device Directive (MDD) 93/42/EEC, Appendix IX.

The manufacturer declares the conformity of this product with the basic requirements according to MDD 93/42/EEC, Appendix VII, and documents this with the CE designation.

English

ITD GmbH
Sportplatzstraße 3
84381 Johanniskirchen
Germany
sales@itd-cart.com
www.itd-cart.com

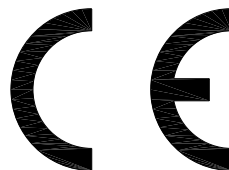Sales, Service and Engineering

ITD GmbH
Grünwalder Weg 13a
82008 Unterhaching
Germany
Telephone: +49 (0)89 61 44 25-0
Telefax: +49 (0)89 61 44 25-20
sales@itd-cart.com
www.itd-cart.com

We work constantly to further develop our products. Please understand that we must reserve the right to make changes to the delivery package in terms of form, equipment and technology at any time.

Reproduction, multiplication or translation, in part or in entirety, is not permitted without the written permission of ITD GmbH!

All rights according to copyright law are expressly reserved for ITD GmbH.

Version © 03/2014

These instructions for usage are valid for the following products:

Part number	Description
CC.02xx.xxx / 03xx.xxx	System components and accessories compact-cart
CC.45xx.xxx / 50xx.xxx / 65xx.xxx	Equipment cart compact-cart , 30 U – 40 U
CC.45xx.xxx / 50xx.xxx / 65xx.xxx	System components and accessories compact-cart
DC.42xx.xxx / 53xx.xxx / 64xx.xxx	Equipment cart duo-cart , 21 U – 30 U
DC.40xx.xxx / 50xx.xxx / 60xx.xxx	System components and accessories duo-cart
EB.45xx.xxx / 50xx.xxx / 65xx.xxx	Equipment cart compact-cart Economy
EC.04xx.xxx	Equipment cart endo-cart , 30 U
GW.52xx.xxx / 53xx.xxx / 54xx.xxx	Equipment cart classic-cart , 21 U – 40 U
GW.62xx.xxx / 63xx.xxx / 64xx.xxx	
GW.01xx.xxx / 02xx.xxx / 03xx.xxx	System components and accessories classic-cart
GW.04xx.xxx / 05xx.xxx / 06xx.xxx	
GW.07xx.xxx / 08xx.xxx	
GW.40xx.xxx / 41xx.xxx / 45xx.xxx	
GW.50xx.xxx / 51xx.xxx / 55xx.xxx	
GW.60xx.xxx / 61xx.xxx / 65xx.xxx	
GW.80xx.xxx / 98xx.xxx / 99xx.xxx	
HA.1xxx.xxx / 2xxx.xxx	System components and accessories flexion-port
HA.45xx.xxx / 5xxx.xxx	
HA.60xx.xxx / 65xx.xxx	
NT.20xx.xxx / 40xx.xxx / 41xx.xxx	System components and accessories vexio-cart
NT.50xx.xxx	Equipment cart vexio-cart , 21 U – 30 U
PT.20xx.xxx / 40xx.xxx / 41xx.xxx	System components and accessories pro-cart
PT.50xx.xxx	Equipment cart pro-cart , 21 U – 35 U
RS.41xx.xxx / 48xx.xxx / 49xx.xxx	Equipment cart uni-cart , 21 U – 35 U
RS.00xx.xxx / 01xx.xxx / 02xx.xxx	System components and accessories uni-cart
RS.4xxx.xxx / 5xxx.xxx	
TS.03xx.xxx / 05xx.xxx	System components and accessories flexion-port

These instructions for usage are valid for the following products:

English

Part number	Description
VS.53xx.xxx / 54xx.xxx / 63xx.xxx	Video rack cart classic-cart , 30 U – 40 U
VS.63xx.xxx	
VT.43xx.xxx	Video cart with isolating transformer pro-cart , 30 U
VT.45xx.xxx	Video rack cart with isolating transformer compact-cart , 40 U
VT.54xx.xxx / 64xx.xxx	Video rack cart with isolating transformer classic-cart , 40 U
VW.54xx.xxx / 64xx.xxx	Video cart classic-cart , 40 U
ZV.94xx.xxx / 96xx.xxx / 97xx.xxx	System components and accessories universal
ZV.98xx.xxx / 99xx.xxx	
KD.7xxx.xxx / 8xxx.xxx / 9xxx.xxx	Customer-specific, mobile equipment carrier, series uni-cart , pro-cart , duo-cart , classic-cart , compact-cart and endo-cart
KU.0xxx.xxx – KU.9xxx.xxx	Customer-specific, mobile equipment carrier, series uni-cart , vexio-cart , pro-cart , duo-cart , classic-cart , compact-cart and endo-cart
OC.0xxx.xxx – OC.9xxx.xxx OC.1 – OC.999	

Content

1	Important information	
1.1	Correct usage	30
1.2	General explanation of the symbols	31
1.3	Safety instructions	33
2	Assembly	
2.1	Completeness	34
2.2	Loading	34
2.3	Loading sequence	34
2.4	Risk through mechanical instability	35
2.5	Castors	35
2.6	Load capacity	35
2.7	Assembly / handling	35
2.8	Additional mounting of system components	36

3	Electrical safety	
3.1	Positioning of the electrical equipment	36
3.2	Power columns (compact-cart , classic-cart , endo-cart) and vertical columns (uni-cart , vexio-cart and pro-cart)	36
3.3	Gases	37
3.4	Equipotential bonding (POAG)	37
3.5	Isolating transformer - leakage current	37
3.6	Earth-Leakage-Guard (ELG)	37
3.7	Cable connection	39
3.8	Combination of equipment	39
3.9	EMC	40
3.10	Excluded from the final electrical inspection of system components and accessories	40
3.11	Minimum safety	41
4	Transportation	
4.1	Safe transportation on castors	41
4.2	Safe transportation by carrying	41
5	Mechanical and electrical height adjustment	42
6	Support arms	
6.1	Wiring	43
6.2	Horizontal swivel	43
6.3	Tilting / rotating equipment	43
6.4	Variable height support arms (flexion-port)	44
6.5	Dismantling and repositioning system components and accessories	44
6.6	Intended use when manoeuvring the mobile equipment carts	44
6.7	Operation of the height adjustable support arms (flexion-port) on the monitor shelf, 2-fold	45
7	Miscellaneous	
7.1	Cleaning and disinfection	45
7.2	Service / Repair	46
7.3	Environmental conditions	46
7.4	Disposal	46
7.5	Spare parts	46
8	Accessories	46
9	Maintenance	47
10	Technical data	
10.1	Load capacity uni-cart	48
10.2	Load capacity vexio-cart	48
10.3	Load capacity pro-cart	48
10.4	Load capacity duo-cart	48
10.5	Load capacity compact-cart	48
10.6	Load capacity classic-cart / endo-cart	48
10.7	Load capacity flexion-port	49

1 Important information

All products from ITD GmbH are manufactured for a long and trouble-free service life. Development, construction, sales and production are certified at ITD GmbH according to ISO 13485.

This is a guarantee for:

- highest quality and a long service life
- easy, safe and ergonomic operation
- functional design
- optimisation for the planned usage

The products satisfy the requirements of the European Medical Device Directive (EMDD) and bear the CE designation.

- Carefully read these instructions for usage from the beginning in order to become familiar with the functions step-by-step.
- Please be sure to address all questions or concerns to the manufacturer.
- The mobile equipment carriers are only intended for the use as described.
- These instructions are to be kept for the service life of the product.

The system configurator is to make the instructions for usage of the overall configuration available to the final customer.

We expressly note here that the system configurator is responsible for the observance of IEC 60601-1 and the EMC norm IEC 60601-1-2 in the valid version!

1.1 Correct usage

The functions of the mobile equipment carriers of ITD GmbH are:

- exclusively to accept medical appliances in accordance with the approved loading data, in compliance with the current applicable edition requirements of IEC 60601-1.
- Connection and distribution of mains voltage from the local electricity outlet as well as from data lines.
- Consolidation of original ITD system components and accessories.

Using the mobile equipment carrier, the medical equipment can be transported inside the building or arbitrarily positioned in the area before and after the application. Therefore, a flexible, economical efficiency of all equipment is possible. In addition, it is easier to clean the floor area.

1.2 General explanation of the symbols

Equipotential bonding: Identifies equipotential bonding terminals on the housing of the isolating transformer; equipotential bonding ensures that resistance between all conductive materials is sufficiently low.

Connection to protective conductor:
Connects conductors, equipment units, conducting parts, main earth terminals and earth


Conductive castors:
Conductive castors are indicated by a flash or a yellow point


Observe the user instructions


Alternating current

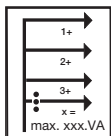
Total output:
Total power that can be supplied by the individual plug in locations must not exceed the total power.

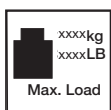Total load rating (basic frame)
max. total load rating (= total of all load ratings of mounted system components)
Please refer to the label for appropriate load.

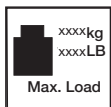Load rating (system components)
Please refer to the label for appropriate load.

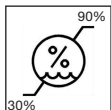Humidity limit

General warning sign:
This symbol is used at the socket strips. The overall rating given on the nameplate must not be exceeded.

Only suitable for the interior

Heavy object:

It is recommended that two people lift the **uni-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** to prevent injury.

Distributed by


Adjustment of the clamping force (tilt and swivel unit)

Set load:

Describes the load range and also the direction of rotation for setting them (3-10 kg and 8-14 kg).


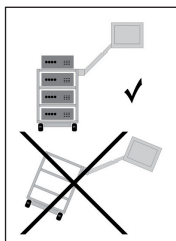
Risk of tilting:

It is imperative to observe the sequence of loading and unloading mobile equipment carts.


Brake:

Information about positioning of the support arm in the event of device acceptance and indicates the direction for release and blocking the locking function.

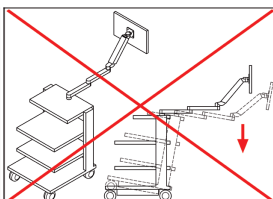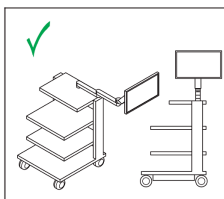Operating position:

This symbol describes the permitted operating positions (right / left) and indicates a prevailing danger of tilting when changing over.

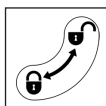Locking function:

This describes the closing / opening direction for components.

Do not push:

It is forbidden to push the equipment cart above the handle because of a risk of tipping over.

1.3 Safety instructions

General

- Only those mobile equipment carriers may be operated whose main voltage equipment has been tested and approved by appropriate, qualified personnel!
- Personnel (hospital and service personnel) working directly or indirectly with a mobile equipment carrier must be instructed!
- Setting adjustments may only be carried out by qualified personnel.
- Repairs may only be carried out by ITD GmbH.

Operation

- With every change of location it must be ensured that no one is injured or objects damaged!

Connections

- In the USA and Canada, use a hospital-grade connection cable when connecting isolating transformers in the 115 V mode and in Japan use the optional Japanese connection cable.
- When connecting to the multiple-socket strip, secure the plug using an optional available plug pull-out lock.
- Only equipment satisfying the requirements of IEC 60601-1 may be connected to the sockets / connecting cables!
- Additional medical equipment with connecting bolts for equipotential bonding are to be connected with the green-yellow cable to the optional equipotential bonding connecting bolts!

- **Caution:** The overall rating given on the nameplate must not be exceeded. Please observe that no further multiple sockets must be connected to an existing multiple socket.

Load capacity

- The total weight of the equipment and the accessories on the mobile equipment carrier may not exceed the permitted payload weight (see load capacity sticker on the base frame).
- The surface load imprinted on the system components may not be exceeded!
- The load shown on the fittings (e.g. infusion stand, jointed arms) may not be exceeded!

Infection protection

- Hygiene regulations are to be observed when cleaning!
- Give only cleaned and disinfected equipment and fittings to a service technician for maintenance and repair work!

Environmental protection

- Dispose of all cleaning and disinfection agent residue in a manner not harmful to the environment!

2 Assembly

2.1 Completeness

Unpack the **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system and check whether all the parts you have ordered have been included in the consignment.

2.2 Loading

Place the **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system on a level, horizontal surface. Place the appliances in the cart taking into consideration the loading sequences. Loading and assembly may only be effected if the cart is separated from current. We recommend the optional accessory (e. g. tension belt) to be applied for the safety of each appliance.


The maximum loading is specified in section 10. If the tilt stability is not adequate, the cart must be clearly labelled in conformity with the 5° labelling specification.

2.3 Loading sequence

Make sure that all the equipment placed on the **uni-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system is stowed securely to prevent slipping, tipping, falling down, etc. by taking appropriate measures (especially also during transport). It is advisable to position all heavy components on the **uni-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system using two people. You should take into account the fact that the centre of gravity changes with the loading.


The cart should be loaded in the following sequence:

- shelves and drawers from the bottom to top
- support arms should be loaded last.

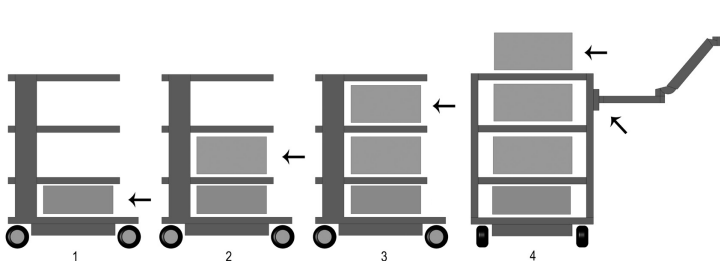The cart should be unloaded in the following sequence:

- loads should be removed from support arms first
- shelves and drawers should be unloaded from top to bottom.

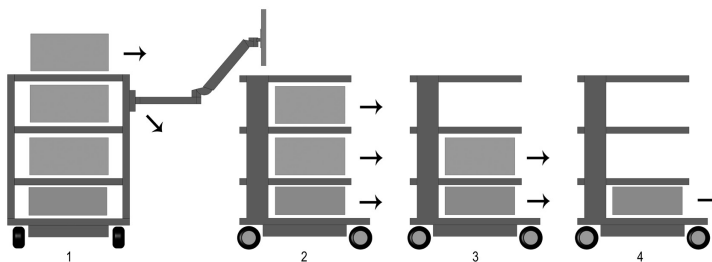
For transportation please refer to section 4 of this manual.

2.4 Risk of mechanical instability

The overall system must be in accordance with IEC 60601-1.

2.5 Castors

The **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system has twin swivel castors with brakes. You should check to see that the wheel locks are working properly before using the **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system. When parking the cart or stopping during transport, all roller brakes on the equipment cart must be pressed (roller locking device). Accordingly, all roller brakes must be released before moving and transport.

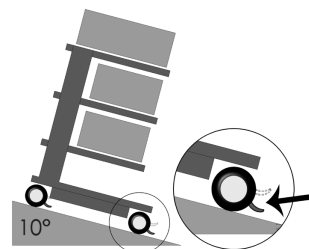The castor must be tested every 12 months for safety as well as the crackproof, firm hub of the castor attachment pin. If these should be lost, please contact your supplier immediately.

2.6 Load capacity

It is not permitted to exceed the load rating of the **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system.

You should observe the maximum load capacity of the **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system (see section 10 of these operating instructions).

2.7 Assembly / Handling


2.7.1 Shelves

Shelves can be removed or installed in another position. Loosen the screws, re-position the shelf and then screw the shelf tightly in position again. Check the safety-earth resistance.

2.7.2 Drawers

Drawer unit blocks (**pro-cart**) are provided with a latch mechanism.

To release the lock on the **pro-cart**, the operating handle on the front cover must be pulled upwards.

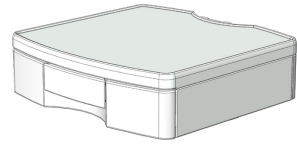The drawers can be lifted out in the extended position. A label strip can be affixed to the front trim (except **pro-cart**). Insert drawer during transport.

English

2.7.3 The isolating transformer

The isolating transformer is mounted in a housing under the base of the basic frame. This is assembled at the factory.

2.8 Additional mounting of system components


Retrospective installation of ITD system components must only be carried out by specialists in accordance with the specifications of the installation instructions supplied. The modified overall system must be rechecked in accordance with the IEC 60601-1.

3 Electrical safety

3.1 Positioning of the electrical equipment

Make sure that electrical equipment on the **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system does not get wet.

Under no condition should you position the products, which may lead to loss of their fluidness, over the electrical equipment or extension lead in which fluidity may be permeate.

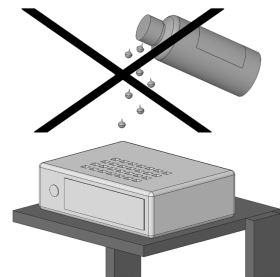3.2 Energy column (**compact-cart**, **classic-cart**, **endo-cart**) and vertical column (**uni-cart**, **vexio-cart**, **pro-cart**)

The socket strips or the cabling system are located in the left or right power columns of **classic-cart**, **compact-cart** or **endo-cart**, with **uni-cart** and **pro-cart** it is in the vertical column. The **vexio-cart**'s socket strips and cabling system are located beneath the base (cabling can be routed through optional cable ducts, which are attached laterally to the vertical column).

Power columns are located on the left and right behind the vertical extrusion and permit optimum accommodation of the power equipment cables available. Do not drill holes in the power columns / vertical columns because power cables may be located inside the columns.

3.3 Gases

Electrical equipment should not be operated in the vicinity of gases, e.g. flammable gas used in anaesthesia or similar gases. The user is responsible for maintaining this requirement and for compliance with IEC 60601-1-2 standard.

3.4 Equipotential bonding (POAG)

Equipotential bonding should be carried out for the **uni-cart**, **pro-cart**, **classic-cart**, **compact-cart** or **endo-cart** carrier system with an isolating transformer. **vexio-cart** offers an optional POAG set without isolating transformer. The equipotential bonding cable should be first connected to the base of the **uni-cart**, **vexio-cart**, **pro-cart**, **classic-cart**, **compact-cart** or **endo-cart** carrier system and then be connected to the equipotential bonding plug in the room. Next, connect the POAG supply cables to the POAG pins of the multiple sockets and the appliance.

3.5 Isolating transformer – leakage current

The **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system is intended to provide a practical and mobile workstation for medical electrical equipment. Total earth leakage currents must not exceed a maximum value of 0.5 mA so that the entire medical electrical system can be compliant with IEC 60601-1. If the sum of the ground leakage currents would exceed the tolerance range, the system must be supplied through a safety transformer.

If an isolating transformer is not mounted, the multiple socket strip / auxiliary socket strip system of the **uni-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system must not be used to connect up equipment that is not in conformity with requirements for earth leakage currents under IEC 60601-1.

If an isolating transformer is mounted, the total power uptake of all the connected units must remain within the rated value for the isolating transformer.


3.6 Earth-leakage-guard (ELG)

The following descriptions are only applicable for the carrier systems in which the item “isolating transformer with insulating monitor” is installed.

3.6.1 Control panel of the earth-leakage-guard (ELG)

At the devices with integrated-earth-leakage guard the evaluation unit is in the transformer case, the separate control panel is placed in one of the shelves. Both components are connected with an interface cord, which lays in one of the profiles.

Control panel of the earth-leakage-guard (optional):

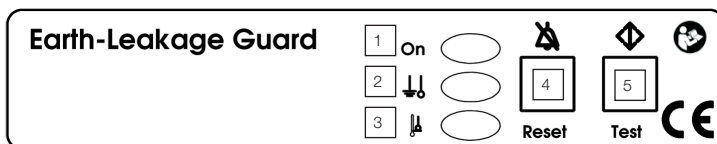- 1 Power on control (green)
- 2 Signals exceeding dielectric resistance (yellow)
- 3 Signals exceeding temperature (yellow)
- 4 Reset-button for fault release
- 5 Test-button starting self-test function manually

3.6.2 Application requirements

The earth-leakage-guard monitors the dielectric resistance for equipment (where safety isolation is necessary) which are connected to the isolation transformers. Simultaneously, the temperature of the transformer is evaluated by a built-in thermal switch which responds when triggered. The evaluation procedure is processor controlled.

3.6.3 Application requirements

Turn the main power switch of the isolating transformer to "ON". Within 5 s the self-test of the ELG is automatically activated.

After the self-test of the ELG is concluded, the green LED (Power on control) glows permanently. During operation this test is repeated (automatically) in a time cycle of 8 hours. This self-test may also be activated manually through the test-button.

The control procedure for the manual self-test is as follows:

- An insulation error is simulated, the yellow LED (signals exceeding dielectric resistance) glows, a warn signal (2.4 kHz) is activated. Both signals stop after approx. 5 s.
- Following this, a temperature error is simulated. The yellow LED (signals exceeding temperature) glows, a warn signal (2.4 kHz pulse tone) is activated. Both signals stop after approx. 5 s.

An error situation is recognized as follows:

- In case of an insulation error, the yellow LED (signals exceeding dielectric resistance) glows permanently, the warn signal with 2.4 kHz is activated permanently. The warn signal can be deactivated through the reset-button for fault release. The LED glows until the error is eliminated.
- If the earth-leakage-guard is turned "OFF" in the meantime and the error is not eliminated, when turning the ELG back "ON" the above described procedure begins anew. With the occurrence of an isolation error the acoustic and optical alarm stays activated (even if the isolation error is corrected on its own) until deactivated through the fault release:

First activation of the fault release: acoustic alarm is deactivated

Second activation of the fault release: optical alarm is deactivated

- In case of a temperature error, the yellow LED (signals exceeding temperature) glows permanently and the warn signal with 2.4 kHz is activated intermittent. The warn signal can be deactivated through the reset-button for fault release. The LED glows until the error is eliminated. If the earth-leakage-guard is turned "OFF" in the meantime and the error is not eliminated, when turning the earth-leakage guard back ON the above described procedure begins anew.
- Should an insulation error and a temperature error occur simultaneously, the dielectric resistance fault has priority (acoustic).

Self test for the Earth-Leakage Guard

In addition to the manual self-test, during operation the earth-leakage-guard repeats the self-test in a time cycle of 8 hours and each time the unit is turned "ON".

From outside, you will not notice the self-test, the time period for this process is approx. 5 seconds. In case of fault, the green LED (power on control) blinks in a frequency of 0.5 Hz and with the same frequency an acoustic alarm is activated. The fault signal can not be deactivated through the reset-button.

3.6.4 Trouble shooting

Do not attempt to repair on your own. Warranty will be invalid if unauthorized repair has been carried out. Due to safety reasons repair and maintenance may only be done by the manufacturer.

Please note:

For further technical data and information, please refer to the separate and comprehensive instructions for use, enclosed with the product, isolating transformer and insulating monitor.

3.7 Cable connection

Users of the **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier system **without an isolating transformer** must ensure that the detachment of the cable plug connection between the cart's multiple power outlet strip and devices can only occur using tools. Our line of accessories includes a selection of covers for multiple power outlet strips.

3.8 Combination of equipment

The following should be observed for combination of equipment on the carrier system:

- Auxiliary equipment connected to analogue and digital interfaces of the equipment must be certified in compliance with the relevant EN specifications (e. g. IEC 60950 for data-processing equipment and IEC 60601-1 for medical electrical equipment).
- Furthermore, all configurations must comply with the valid version of the standard IEC 60601-1. Anyone connecting additional equipment to the signal input or signal output is a system configurer and is therefore responsible for ensuring compliance with the valid version of the standard IEC 60601-1.

If you have any questions, you should contact your local dealer or technical services.

Please note:


Make sure that this is also the case for adaptation of equipment in the power supply circuit (e. g. multiple socket strip).

3.9 EMC


Electromagnetic compatibility of the medical electrical equipment positioned on the cart has to be checked by the overall system configurer. If you use a different equipment combination, you should check electromagnetic compatibility between individual items of equipment before using in a medical setting.

The carrier system cannot be used in an Magnetic Resonance environment due to the presence of ferro-magnetic materials.

ITD GmbH excludes any liability in this respect.

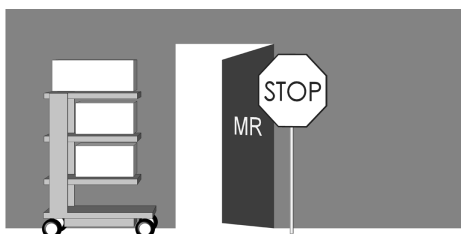Using the carrier system with an isolating transformer within an environment subject to explosive hazard is not permissible.

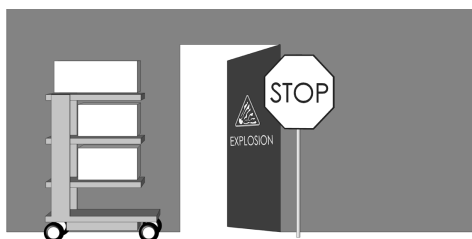3.10 Excluded from the final electrical inspection of system components and accessories

ITD GmbH exclude the following system components and accessories from the final electrical inspection:

- Multiple sockets that are not wired in the assembly
- ME cables and appliance cables included
- POAG plates and cables included

- Non-electrified equipment carts and support systems
- Height adjustments and attachment parts to height adjustments
- Handles, mouse pads, drawers, drawer bodies and attachments (bottle mounting brackets, baskets, camera mounting brackets, infusion tripods, ...)
- Keyboard extensions and extendable shelves
- Computer holder
- Conductive castors
- Support arms installed and monitor mounting brackets

„Secondary power circuit with insulation monitors“ are only excluded from the dielectric strength inspection!

3.11 Minimum safety

ITD GmbH is not aware of any item of equipment or accessory that reduces the minimum safety of the system. Only equipment not presenting a hazard may be used.
If necessary, this should be clarified by means of a risk analysis (ISO 14971).

4 Transportation

4.1 Safe transportation on castors

Please ensure that before the mobile use of the equipment cart:

- All equipment / products in transit are secured against detachment.
- All swivel arms must be swivelled in and secured.
- The feeder from the local outlet must be disconnected.
- The brakes of the castors are released.

The equipment cart should only be moved at a maximum walking pace of 6 km/h. The velocity is significantly reduced when knocked down by bedplates, cables as well as hoses. It is safer to deliver by driving on ramps having a max. gradient of 10° so that the equipment cart can be brought to a standstill at any time.

Although all precautionary measures have been observed to guarantee maximum stability for this product, it is important to pay attention to unevenness in floors, lift door frames, cables, etc. in order to prevent accidents.

Fundamentally, the requirements of the IEC 60601-1 apply.

4.2 Safe transportation by carrying

Handles are not provided for lifting the mobile equipment carts, but exclusively have a function for pushing. Lifting and carrying the mobile device carrier must only be carried out by two persons using the extensions of the base.

Fundamentally, the requirements of the IEC 60601-1 apply.

5 Mechanical and electrical height adjustment

Specific safety regulations must be observed for the „mechanical height adjustment“ using the „gas pressure“, as well as for the electro-mechanical height adjustment using „Linear drive“, in accordance with IEC 60601-1 „Mechanical Risk in Conjunction with Moving Parts“. Hereby:

- Take into consideration and adhere to the permitted distance between moving parts, in accordance with IEC 60601-1 in Table 20 (ISO 13857:2008).
- Products with height adjustment are manufactured and supplied ex works conforming to the standard, under consideration of the permitted safety distances. Due to the equipping of or replacement with ME devices and / or component, these distances change. This can result in a mechanical risk. The respective person who configures the system is responsible for adhering to the minimum distances required.
- The overall weight of the devices and accessories installed must not exceed the specified maximum overall payload of the height adjustment. Overloading results in damage to the height adjustment and loss of the warranty.
- Stored energy is released with the mechanical height adjustment using gas pressure. Thereby, for unloaded systems sudden, unbraked activation of the height adjustment can result in injuries and damage.
 - o In order to prevent injuries and damage, before installation and removal of the devices, place the height adjustment at the uppermost position („energy-free“).
 - o The height adjustable supporting arm system „**flexion-port**“ must also be fixed and secured at the uppermost position („energy-free“) with the aid of the clamping lever (refer to the separate user instructions for the „**flexion-port**“, as well as the risk notes on the supporting arm system).
- Unintended activation of the electro-mechanical height adjustment using the manual button can also result in injuries and damage.
 - o In order to prevent injuries and damage, before installation and removal of the devices, disconnect the height adjustment from the power supply.
 - o Servicing and maintenance tasks in the „interior area“ of the height adjustment, i.e. in the covered area within the supporting column not accessible from outside, must only be carried out by specialists.
 - o **Caution:** If operating the height adjustment using a remote control, make sure that no persons are in the hazard area.

6 Support arms

6.1 Wiring


Please observe the following instructions:

- So that there is no damage to the cable or failure of the device when swivelling, the cable must be sufficiently dimensioned.
- Possible sagging cables must under no circumstances be used as a handle.
- Please make sure that the enclosed assembly material is correctly applied, in accordance with the assembly instructions.
- When swivelling the arms, pay attention to possible cable loops present.

6.2 Horizontal swivel


Make sure that the swivel range of the system components matches the dimensions of the equipment and the ambient conditions in the working environment.

When system components and attached equipment are swivelled horizontally, they must not collide with other equipment, with other system components or with the wall. Any collision may result in damage to equipment and injuries to persons.

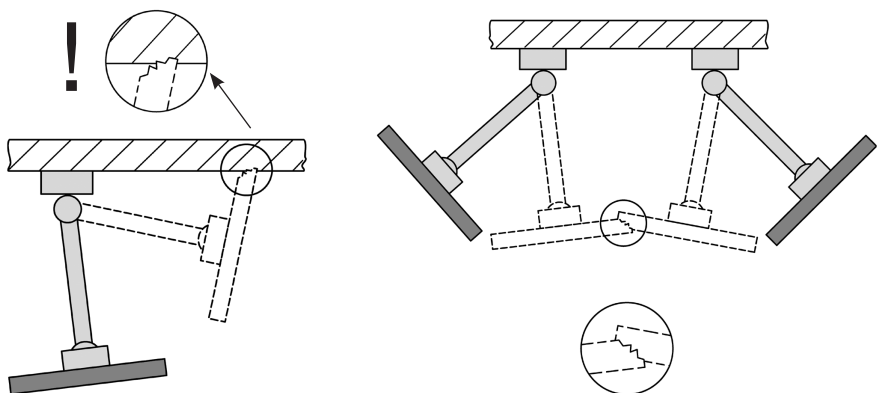
6.3 Tilting / rotating equipment

If system components are tiltable or rotatable, it is important to ascertain whether the clamping force is appropriate for the unit being fixed in place. If the force applied is incorrectly adjusted, the equipment is at risk of tilting over. Adjustment must therefore be carried out to ensure that the unit can be slightly tilted or rotated, while the unit remains stable in any desired position.

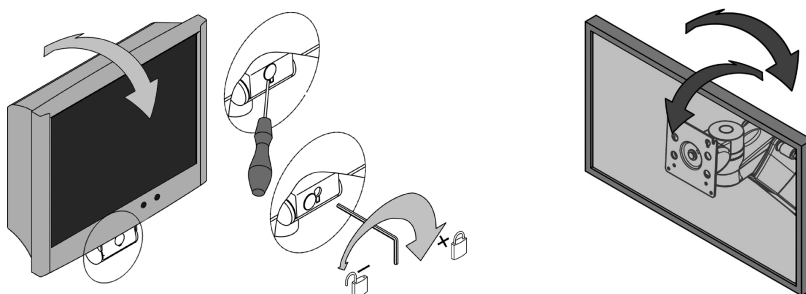6.4 Variable height support arms (**flexion-port**)

When loading the system components that can be height adjusted, unconditionally observe the minimum and maximum permitted total weight. Also, due to safety reasons, please make sure that the space below the height-adjustable support arm (**flexion-port**) remains clear. In order to be able to set the support arm to the load, it must be placed into the horizontal position.

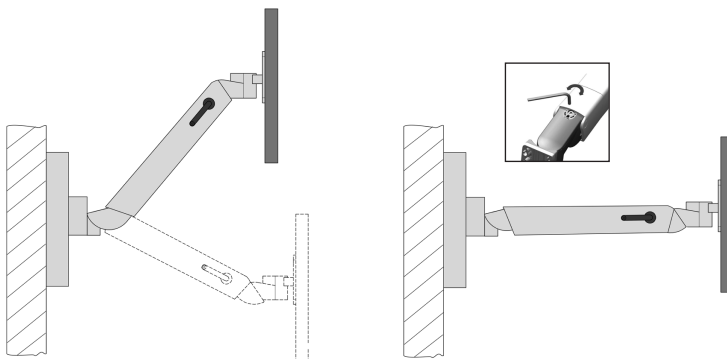6.5 Dismantling and repositioning system components and accessories

When system components and accessories are dismantled or repositioned, it is important to remove any units mounted on these elements before any changes are made. If this concerns the disassembly / assembly of the height-adjustable support arms **flexion-port**, these must first be placed in the uppermost position and the clamping (brake) must be determined (refer to decal).

6.6 Intended use when manoeuvring the mobile equipment carts

When manoeuvring the mobile equipment carts it is imperative to make sure that the support arms attached (swivelling, height adjustable) are placed into the respective park position and, as necessary, locked. Otherwise, the stability cannot be guaranteed (refer to Point 2.4).

6.7 Operation of the height adjustable, support arm system (**flexion-port**) on the monitor shelf, 2-fold

When operating the height adjustable support arm system (**flexion-port**) on the monitor shelf, 2-fold, always make sure that the swivelling arm that provides the changeover, is always placed parallel to the front side of the cart (engaged position). Thereby, please observe that only two operating positions of the support arm (right / left) are permitted (refer to „Operating position“ in Chapter „1.2 General symbol explanation“, Page 31).

By a changeover of the operating position from left to right, or vice versa, the locking elements must be pushed downwards and the height adjustable support arm swivelled to the other side. Hereby, for **flexion-port** with swivelling arm make sure that it is folded when changing over. Otherwise, the tilting stability cannot be guaranteed (refer to Point 2.4).

7 Miscellaneous

7.1 Cleaning and disinfection

Caution: Disconnect from power before cleaning and disinfection!

The equipment cart **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** must be cleaned using commercial all-purpose cleaning agents (neutral cleaning agents).

For disinfecting, commercial disinfectants approved for disinfecting surfaces or wipe disinfection can be used. The disinfection agents must be used solely as disinfection for wiping, in accordance with the manufacturer's specification.

For example, ITD have carried out tests using the following disinfectants:

Product	Manufacturer
Bacillol plus	Bode
Cleanisept Wipes	Dr. Schumacher
Mikrobac Tissues	Bode
Mikrozid Sensitive Wipes	Schülke
Terralin Protect	Schülke
Incidin PLUS	Ecolab
Incidin Foam	Ecolab

If complete disinfection is required, assemblies can be disassembled by a specialist and wipe disinfected when disassembled.

7.2 Service / Repair

The cart should be always be cleaned and disinfected with a suitable cleaning agent before any service operations are undertaken, and before the cart is returned for purposes of repair. Repairs to the cart should only be effected by professional personnel. We recommend consulting ITD GmbH on all matters relating to service activities.

ITD GmbH
Grünwalder Weg 13a
82008 Unterhaching, Germany
Tel: +49 (0)89 61 44 25-0, Fax: +49 (0)89 61 44 25-20
sales@itd-cart.com, www.itd-cart.com

7.3 Environmental conditions

The **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** carrier systems are designed for standard operation in hospitals and medical practices.

Ambient temperature: 10° C to 40° C
Air humidity: 30 % to 75 %
Air pressure: 700 hPa to 1060 hPa
Protection class: IP20

7.4 Disposal

Separate Collection for Electrical and Electronic Equipment in compliance with Waste Electrical and Electronic Equipment Directive WEEE (registration number for Germany: DE35464575). All electrical and electronic equipment provided with systems released after 13 August 2005 is marked with a Separate Collection for Electrical and Electronic Equipment symbol, indicating that this equipment must undergo separate collection for disposal, in countries where EU directive 2002/96/EC is in effect.

7.5 Spare parts

Only spare parts authorized by ITD may be used. A sticker with an order number is attached to the base of your cart. All order numbers and the associated spare parts are archived at ITD GmbH. Spare parts may be obtained from the following address:

ITD GmbH
Grünwalder Weg 13a
82008 Unterhaching, Germany
Tel: +49 (0)89 61 44 25-0, Fax: +49 (0)89 61 44 25-20
sales@itd-cart.com, www.itd-cart.com

8 Accessories

A comprehensive range of accessories is provided in our catalogues or under www.itd-cart.com (information for dealers).

9 Maintenance

The **uni-cart**, **vexio-cart**, **pro-cart**, **duo-cart**, **compact-cart**, **classic-cart** or **endo-cart** has been developed and constructed for many years of trouble-free use.

Check the functional capability of the following parts every 12 months in order to guarantee safety.

Monitor shelf:

- Swivels and tilts easily without too much play.

Shelves:

- Check whether the mounting screws have been tightened and whether the shelf is stable and flat.

Castors:

- Check that the castors run freely and the brakes operate correctly.
- Check that the 4 bolts holding the castors to the bottom side of the base and the castors themselves are seated firmly in their holders.
- The running surfaces of castors must be free from contamination, so that they can function correctly.

Socket strips:

- Check the main cable for damage and firm seating.

Auxiliary sockets:

- Check the cable for damage and firm seating.

Variable height support arms (**flexion-port**):

- The height adjustment functions freely, the raising force is adjusted to the weight of the device.

Support arms:

- Rotating and tilting functions smoothly, without excess play.

Isolating transformers:

- Safety-relevant checks of isolating transformers.

Serial number:

- Compare the serial number of the cart with the data of the equipment log book.

Fuses:

- Check whether the correct fuses have been installed.

If you encounter any problems during these checks you should contact your supplier immediately.

10 Technical data

10.1 Load capacity **uni-cart**

- Basic frame, total payload from 50 kg / 110 lbs
- Shelf 10 kg / 22 lbs
- Rack shelf 20 kg / 44 lbs
- Drawer 3 kg / 6.6 lbs (lockable: 10 kg / 22 lbs)
- Monitor holder 14 kg / 30.8 lbs

10.2 Load capacity **vexio-cart**

- Basic frame, total payload 65 kg / 143 lbs
- Shelf 15 kg / 33 lbs
- Drawer unit 3 kg / 6.6 lbs
- Monitor holder 14 kg / 30.8 lbs

10.3 Load capacity **pro-cart**

- Basic frame, total payload 80 kg / 176 lbs
- Shelf 20 kg / 44 lbs
- Drawer unit 15 kg + 3 kg / 33 lbs + 6.6 lbs
- Monitor holder 14 kg / 30.8 lbs

10.4 Load capacity **duo-cart**

- Basic frame, total payload 80 kg / 176 lbs
- Shelf 50 kg / 110 lbs (pull-out: 20 kg / 44 lbs)
- Drawer unit 3 kg / 6.6 lbs (lockable: 20 kg / 44 lbs)

10.5 Load capacity **compact-cart**

- Basic frame, total payload 180 kg / 396 lbs
- Shelf 50 kg / 110 lbs
- Drawer unit 3 kg / 6.6 lbs
- Monitor shelf max. 35 kg / 77 lbs (depends on type)

10.6 Load capacity **classic-cart / endo-cart**

- Basic frame, total payload 150 kg / 330 lbs
- Shelf 50 kg / 110 lbs (pull-out: 20 kg / 44 lbs)
- Drawer unit 3 kg / 6.6 lbs
- Monitor shelf max. 35 kg / 77 lbs (depends on type)

10.7 Load capacity **flexion-port** (variable height support arms)

- **flexion-port** 3-10 kg / 6.6-22 lbs or
8-14 kg / 17.6-30.8 lbs
- Tilt and swivel unit /
rotation unit 14 kg / 30.8 lbs
- Post ("Down-Post") 10 kg / 22 lbs

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ITD GmbH si assume la responsabilità dell'esattezza del contenuto.